

# Magnetic Anisotropy: Measurements with a Vector Vibrating Sample Magnetometer

B. C. Dodrill, J. R. Lindemuth, and J. K. Krause

*Lake Shore Cryotronics, Inc., 575 McCorkle Blvd., Westerville, Ohio, USA*

An automated vector vibrating sample magnetometer (VVSM) is described which, in addition to standard hysteresis loop measurements, provides for measurement of the angular dependence of the vector components of the total magnetization. This is of particular utility in investigations of rotational hysteresis and anisotropy in magnetic materials.

The VVSM affords the opportunity to determine anisotropy constants from torque curves. Two sets of sensing coils are employed to measure the components of  $\mathbf{M}$ , from which  $M_{\parallel}$  the component parallel to the applied field  $H$ , and  $M_{\perp}$  the component perpendicular to  $H$  are deduced. The torque density exerted on the sample by the applied field is  $\tau = \mu_0 \mathbf{M} \times \mathbf{H} = -\mu_0 M_{\perp} H \hat{\mathbf{K}}$ . Data are presented for metal particle and  $\text{CrO}_2$  tape media with in-plane anisotropy, and compared with results obtained from a torque magnetometer.

## Introduction

Anisotropy is a basic property of magnetic materials<sup>1,2,3</sup>. The magnetization tends to lie in certain preferred directions, and the energy of the ferromagnet includes an anisotropy term  $E_a(\theta, \phi)$  where the angular coordinates  $\theta, \phi$  define the direction of magnetization. The anisotropy can be intrinsic, related to atomic scale interactions which define easy directions in the crystal (magnetocrystalline anisotropy), or it can be related to the energy of the sample in its own demagnetizing field (shape anisotropy). Magnetic anisotropy is a very important parameter in connection with the characterization of materials utilized in technical applications, particularly magnetic recording media.

The usual method for anisotropy determinations is torque magnetometry, which directly measures the macroscopic torque exerted by an applied field on the magnetization of a sample that is not precisely aligned with the field. Although torque magnetometers enjoy the advantage of high sensitivity and accuracy, they are unsuited to measure other magnetic parameters of interest, e.g., hysteresis loops, remanence curves, etc. Measuring magnetization curves for different orientations of a sample with a normal Vibrating Sample Magnetometer (which measures only  $M_{\parallel}$ , the parallel component of the magnetization) is an indirect method of studying the anisotropy field. Here we present a vector magnetometer which is capable of measuring the magnetization component perpendicular to the applied field. This component is directly related to the macroscopic torque<sup>4,5,6,7,8</sup>, and hence the VVSM provides information that is essentially identical to that provided by a torque magnetometer. In addition, the VVSM offers all the advantages of a normal VSM for measuring magnetization, susceptibility, hysteresis, and remanence<sup>9</sup>.

## Theory: Anisotropy Constants from Vector Magnetization Data

The anisotropy is obtained in a manner very similar to torque magnetometry, the difference being that instead of measuring torque directly, which is the product  $\mathbf{H}M_{\perp}$ , one obtains the transverse magnetization  $M_{\perp}$ , and the longitudinal magnetization  $M_{\parallel}$  from the components  $M_x$  and  $M_y$  measured by the vector magnetometer coils. These experimental quantities are related to the torque and anisotropy constants. The formalism is as follows.

The energy of the system is given by the sum of anisotropy and magnetic potential energies,

$$E_T = E_a + E_M$$

In the case of uniaxial symmetry we have,

$$E_a = K_0 + K_1 \sin^2 \theta + K_2 \sin^4 \theta + \dots$$

Where  $\theta$  is the angle between  $\mathbf{M}$  and the easy direction of magnetization of the sample, and  $K_0$  is constant and independent of angle. The potential energy is given by  $-\mu_0 \mathbf{M} \cdot \mathbf{H}$ , hence

$$E_p = -\mu_0 MH \cos(\psi - \theta)$$

where  $\psi$  is the angle between  $\mathbf{H}$  and the easy direction of magnetization (see figure 1). If  $K_2 \ll K_1$

$$E_T = K_0 + K_1 \sin^2 \theta - \mu_0 MH \cos(\psi - \theta) \quad (1)$$

At equilibrium the total energy is minimized, which requires that  $dE_T/d\theta = 0$ , hence

$$K_1 \sin(2\theta) = \mu_0 MH \sin(\psi - \theta). \quad (2)$$

The expression on the left is the torque  $\tau = -dE_a/d\theta$  exerted on the magnetization  $\mathbf{M}$  by the crystal, and the expression on the right is the torque exerted on  $\mathbf{M}$  by the applied magnetic field. Since  $M \sin(\psi - \theta) = M_{\perp}$ , the transverse magnetization,

$$K_1 \sin(2\theta) = \mu_0 M_{\perp} H.$$

When  $M_{\perp}$  is a maximum,  $\sin(2\theta) = 1$ . Hence,

$$K_1 = \mu_0 M_{\perp}^{\max} H. \quad (3)$$


Measurement of  $M_{\perp}$  as a function of angle  $\psi$  therefore allows determination of the anisotropy constant  $K_1$  and torque  $\tau$ . In the VSM where a resistive electromagnet is employed for field generation, the angle  $\psi$  is varied by rotating the sample relative to the magnetic field that is fixed in direction.

An alternative method of extracting a value for  $K_1$  and also a value for  $K_2$  is by fitting the torque curve to a Fourier series. All data in the torque curve is then used to calculate the anisotropy, not just the peak value.

Keeping terms in the energy to order  $\sin^4 \theta$ , the torque is:

$$\tau(\theta) = (K_1 + K_2) \sin(2\theta) - 2K_2 \sin(4\theta).$$

$K_1$  and  $K_2$  can both be extracted by least-squares fitting the torque curve to this function. If the magnitude of the Fourier components for  $2\theta$  and  $4\theta$  are  $\tau_{2\theta}$  and  $\tau_{4\theta}$  respectively, then  $K_2 = -\tau_{4\theta}/2$ ,  $K_1 = \tau_{2\theta} - K_2$ . If we only keep terms to order  $\sin^2 \theta$  in the energy, then  $K_1 = \tau_{2\theta}$ .

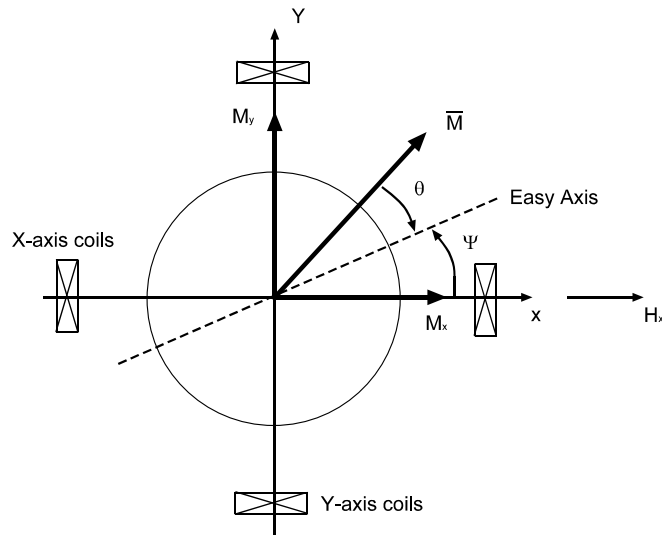

**Figure 1:** Schematic top view of the EM vector VSM


## Experimental

A schematic illustration of an electromagnet vector VSM (EM-VVSM) is shown in figure 1 with the different angles defining the directions of magnetization and applied field. One pair of sensing coils are parallel to the applied field and sense the magnetization component longitudinal to the field  $M_{||} = M_x$ . A second set of coils is mounted at right angles to the applied field and senses the magnetization component transverse to the field  $M_{\perp} = M_y$ . Hence the VVSM is a *biaxial* VSM and may be used to measure anisotropy in the  $xy$  plane. Although not considered in the present discussion, note that a third set of coils to measure  $M_z$  may also be utilized. In the case where the magnetization is a true 3-dimensional vector, e.g., perpendicular recording media, obliquely deposited tapes with anisotropy tilted out of the film plane, etc, the *tri*-axial VSM provides for similar anisotropy determinations. Sample rotation,  $\psi$ , in the  $xy$ -plane is achieved via a computer-controlled motor attached to the VSM head.  $M_x$  and  $M_y$  are recorded simultaneously as a function of  $(H_x, \psi)$  under full computer control in the EM-VVSM.


## Results

Data were recorded for a iron metal-particle (MP) tape sample, and a  $\text{CrO}_2$  tape sample. Tape areas were  $0.317 \text{ cm}^2$ , and tape thickness were  $0.69 \text{ }\mu\text{m}$  and  $0.55 \text{ }\mu\text{m}$ , respectively. Both  $M_x$  and  $M_y$  were recorded as a function of magnetic field  $H_x$ , and as a function of angle  $\psi$  in the EM-VVSM. The  $\text{CrO}_2$  tape was also measured using a conventional torque magnetometer.


Figures 2 and 3 illustrate results of the EM-VVSM measurements and show the response of both the  $x$  and  $y$  axis coils as a function of magnetic field for various angles  $\psi$  of the applied field relative to an initial arbitrary orientation  $\psi_0$  with respect to the sample's easy axis for the iron-particle and  $\text{CrO}_2$  tapes, respectively. The  $M_x(H_x)$  response is the classical hysteresis loop normally measured in a single-axis VSM.

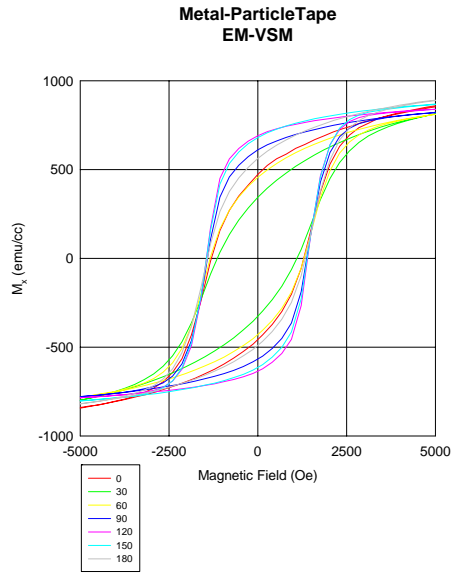


**Figure 2a:**  $M_x(H_x)$  for a MP tape sample

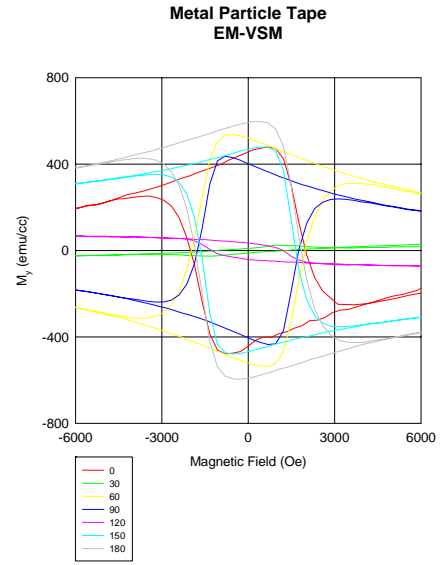

**Figure 2b:**  $M_y(H_x)$  for a MP tape sample

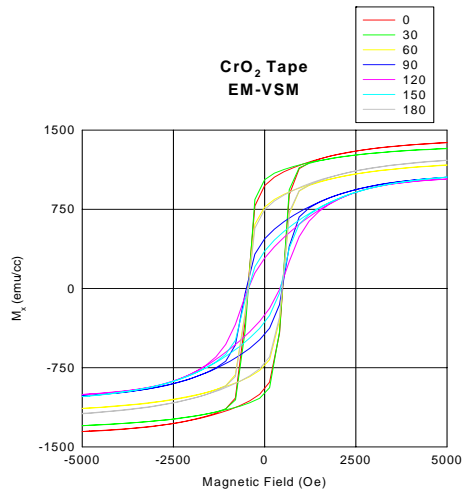


**Figure 3a:**  $M_x(H_x)$  for a CrO<sub>2</sub> tape sample

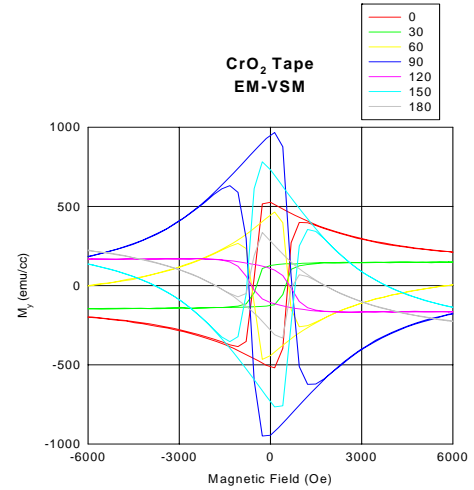


**Figure 3b:**  $M_y(H_x)$  for a CrO<sub>2</sub> tape sample

Shown in figure 4 are the x and y moments recorded at a field strength of 5 kOe as a function of angle from 0° to 360° as recorded on the EM-VVSM for an MP tape sample. The  $m_y$  value goes to zero when the easy axis of magnetization is aligned with the field direction. At remanence,  $m_y = 0$  and  $m_x = \text{maximum}$  when the easy axis is aligned with the x-axis coils, or alternatively  $m_y = \text{maximum}$  and  $m_x = 0$  when the easy axis is aligned with the y-coils. Thus, automated easy axis determinations are possible employing the vector VSM technique.

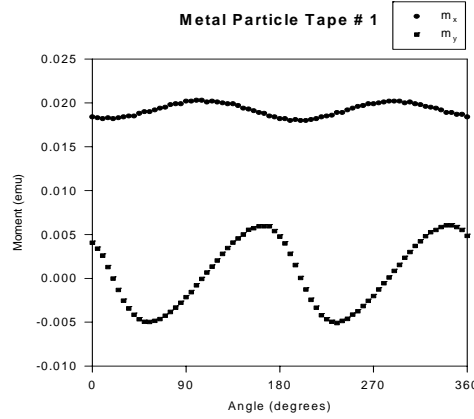

**Figure 4:**  $m_x(\psi_0 + \psi)$  and  $m_y(\psi_0 + \psi)$  at 5 kOe for a MP tape sample

Applying equation 3,  $K_1 = M_y H = m_y H / V$  at the maximum in  $m_y$ . Alternatively, the torque (in dyne-cm)  $\tau = \mathbf{m} \times \mathbf{H} = m_y H_x \mathbf{k}$ , and hence  $\tau(\theta)$  may be similarly derived. From the data shown in figure 4, the anisotropy constant for the ME tape is:  $K_1 = 1.27 \times 10^6$  ergs  $\text{cm}^{-3}$ . From similar data recorded for the  $\text{CrO}_2$  sample,  $K_1 = 0.91 \times 10^6$  ergs  $\text{cm}^{-3}$ . Figure 5 shows the angular dependence of the torque for both the MP and  $\text{CrO}_2$  tapes.

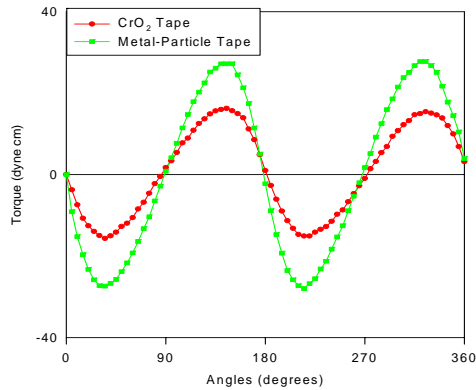


**Figure 5:**  $\tau(\theta)$  (dyne-cm) for both MP and  $\text{CrO}_2$  tapes as determined from the EM-VVSM data

To check the validity of the anisotropy derived from the EM-VVSM data, the CrO<sub>2</sub> tape sample was also measured in a normal torque magnetometer. Figure 6 shows the resultant torque density (dyne cm<sup>-2</sup>) curves for the EM and torque magnetometer. The anisotropy constants determined from these data are:

	Torque Magnetometer	EM-VVSM
<b>K<sub>1</sub> (ergs cm<sup>-3</sup>)</b>	0.96 × 10 <sup>6</sup>	0.91 × 10 <sup>6</sup>

These data agree to the few percent level, thus verifying the vector VSM anisotropy determination. In order to obtain reliable values of the anisotropy constants, especially K<sub>2</sub>, it is important to use an applied field that is sufficient to saturate the sample. A good K<sub>2</sub> value would require at least double the 5.5 kOe fields used in these measurements.

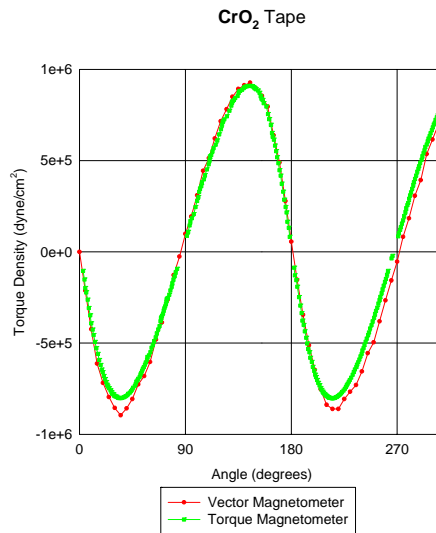

**Figure 6:** Torque density vs. angle for CrO<sub>2</sub> tape as determined from the EM vector VSM and a torque magnetometer

## Summary

This paper has discussed the use of a vectorial measurement employing a bi-axial VSM for evaluating the anisotropic magnetic behavior of metal particle and chromium oxide tape samples exhibiting uniaxial in-plane anisotropy due to alignment of the acicular particles. In addition to measuring hysteresis loops, the VVSM may also be used to measure anisotropy and torque curves. It has been shown that the VVSM yields results in good agreement with curves measured using a classical torque magnetometer. Although the discussion here was limited to a 2-dimensional treatment, extension to 3-dimensional anisotropy, for example in obliquely deposited films with anisotropy tilted out of the film plane and perpendicular magnetic media, is possible employing a tri-axial VSM.

## References

1. R. F. Pearson, in Experimental Magnetism, Vol. 1 (G. M. Kalvius and R. S. Tebble, editors) Wiley, Chichester 1979, Ch. 3.
2. Ferromagnetism, Richard M. Bozorth, IEEE Press, NY 1951
3. Introduction To Magnetic Materials, B.D. Cullity, Addison-Wesley, MA, 1972
4. Magnetometer For Anisotropy Measurement Using Perpendicular Magnetization, E. Joven, A. Del Moral, and J. I. Arnaudas, Journal Of Mag. And Mag. Mat, 83, pg. 548, 1990
5. Analysis of Perpendicular Recording Media Using Biaxial VSM, K. Ouchi and S. Iwasaki, IEEE Trans. On Mag., 24 (6) pg 3009, 1988
6. Measurement of Magnetization on an Obliquely Deposited Film by Using 3-dimensional Vector VSM, N. Matasubara and F. Sai, IEEE Trans. On Mag., 27 (6), pg. 4748, 1991
7. Fourfold Anisotropy of An Electrodeposited Co/Cu Compositionally Modulated Alloy, R. D. McMichael, U. Atzmony, C. Beauchamp, L. H. Bennett, L. J. Swartzendruber, D. S. Lashmore, and L. T. Romankiu, Journal Of Mag. And Mag. Mat., 113, pg. 149, 1992
8. R. Josephs, D. S. Crompton, and C. S. Kraft, Invited Talk, Intermag Conference, 1987
9. Versatile and Sensitive Vibrating Sample Magnetometer, S. Foner, Rev. Sci. Instr, 30, pg. 548, 1959