




## RGC Series recirculating gas coolers

**RGC Series** convert a liquid helium cryostat to operate without liquid helium, for cryogen-free operation. Compatible with Lake Shore ST and STVP cryostats.

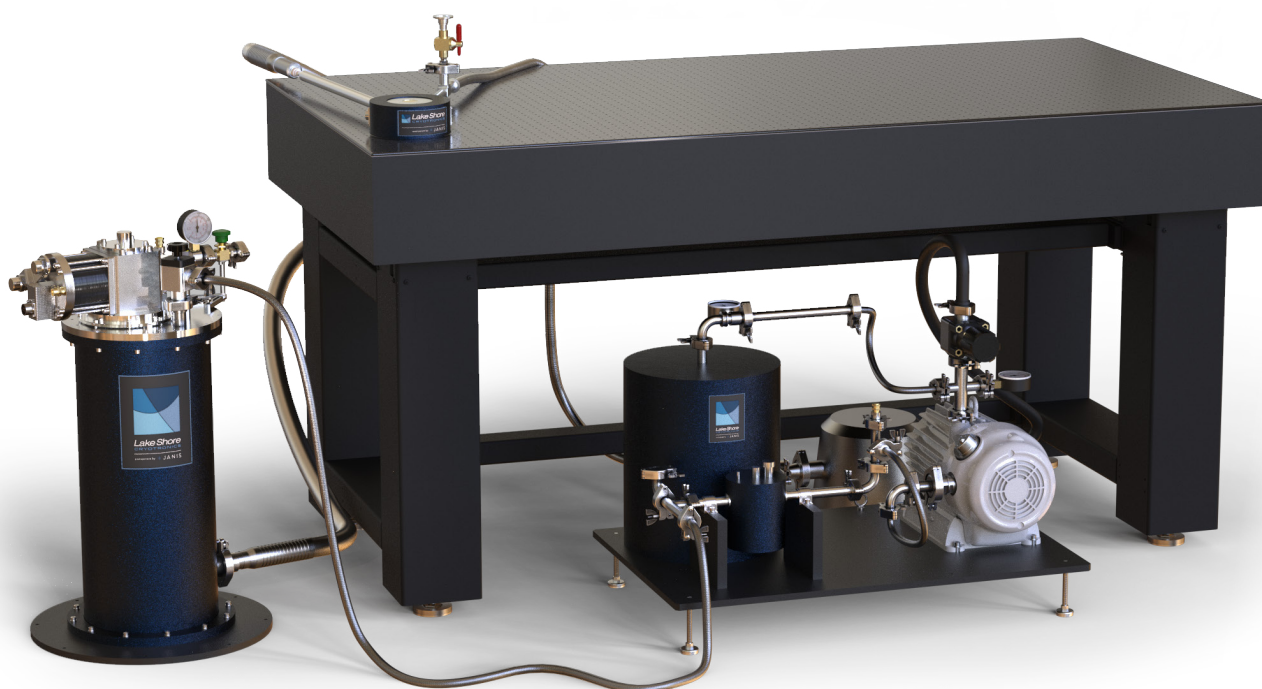

## RGC Series recirculating gas coolers

RGC Series recirculating gas coolers run helium in a closed loop, making a continuous flow cryostat cryogen-free. Helium gas is cooled and liquefied by the RGC's cryocooler, and travels to the cryostat through a flexible vacuum-insulated transfer line. LHe cools the sample. The RGC captures the evaporated gas through the transfer line and reliquefies it, continuously recirculating the helium.

Samples can be changed without warming up the RGC, allowing fast turnaround times. When paired with a Lake Shore ST-500 Series cryostat, the combination is an ultra-stable cryogenic microscopy platform, cooling samples and devices to below 4 K without the use of liquid helium.

The RGC is compatible with Lake Shore ST and STVP cryostats, and can be used with some LHe cryostats from other vendors as well.

### Key features

---

Maintain benefits of a LHe cryostat without the LHe

---

Ideal for low vibration

---

Excellent thermal performance

---

Fast sample change without warming up the RGC

---

## RGC Series

### Featured components

Two-stage Gifford-McMahon (G-M) cold head with compatible compressor (pulse tube cooler optional)


Vacuum shroud containing internal helium lines with special heat exchangers for cooling the circulating helium gas, radiation shield, and two diagnostic standard curve silicon diode temperature sensors

Gas handling system with hermetically sealed oil-free gas circulating pump, storage volume, and interconnecting pumping lines and valves


Vibration isolation for gas circulation line

Helium transfer line to deliver helium to cryostat

Base temperature of ST-500 microscopy cryostat cooled by RGC4 system

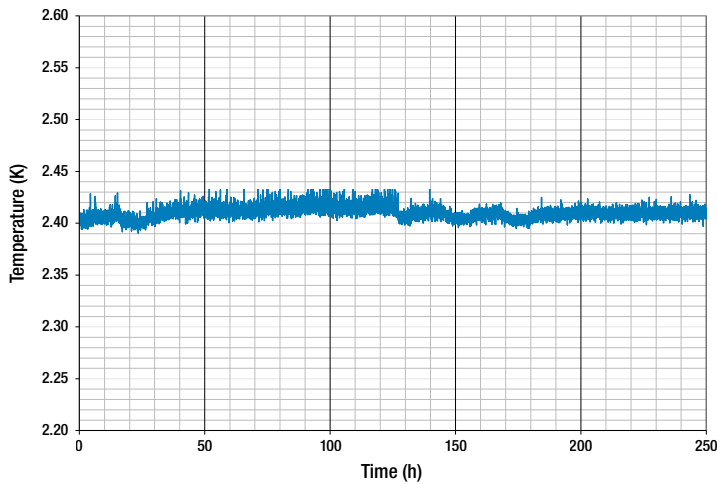

Vibration measured on a standard ST-500 cryostat cooled by an RGC system

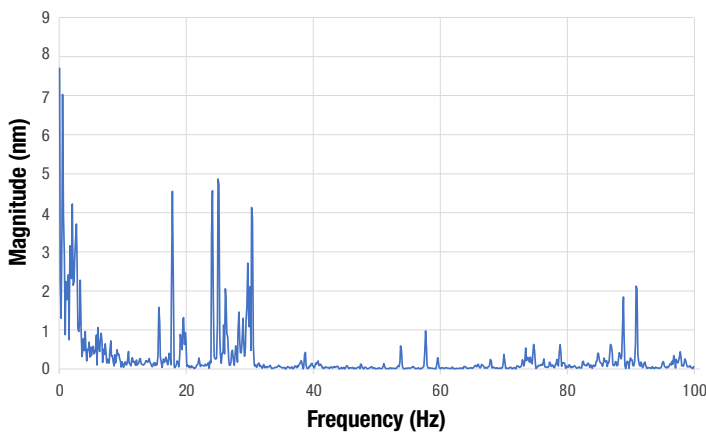

## Selections

### Cryocooler models

RGC4-10: 1 W

RGC4-15: 1.5 W

RGC4-20: 2 W

### Compressor

Water-cooled

Air-cooled

### Gas handling system pump

3.6 cfm scroll pump (standard)

10.6 cfm scroll pump

### Custom configurations

Pulse tube cryocooler instead of G-M cryocooler

Elbow on transfer line

## Compatible cryostats

Compatible with Lake Shore ST and STVP cryostats, and can be used with some LHe cryostats from other vendors as well.

## Specifications

| Sample temperatures <sup>1</sup> | ST-100 and ST-300 | ST-400                 | ST-500                 | STVP  | Probe station      |
|----------------------------------|-------------------|------------------------|------------------------|-------|--------------------|
| RGC4-10                          | <4.3 K            | <4.0 K (120 mW at 5 K) | <4.2 K (100 mW at 5 K) | <10 K | Consult Lake Shore |
| RGC4-15                          | <3.5 K            | <2.9 K (220 mW at 5 K) | <3.8 K (210 mW at 5 K) | <8 K  |                    |
| RGC4-20                          | <3.3 K            | <2.6 K (280 mW at 5 K) | <3.5 K (250 mW at 5 K) | <7 K  |                    |

### Size

|   |   |
|---|---|
| Gas handling system dimensions                                | 533.4 mm × 838.2 mm × ≈711.2 mm (21 in × 33 in × ≈28 in)  |
| RGC cryocooler dimensions                                     | 457.2 mm (18 in) base diameter × 787.4 mm (31 in) tall  |
| Gas handling system weight (approximate)                      | 56.7 kg (125 lb)  |
| RGC cryocooler weight (approximate)                           | 36.3 kg (80 lb)   |
| Shipping dimensions for three RGC system crates (approximate) | Crate 1: 109 cm × 84 cm × 147 cm (43 in × 33 in × 58 in);<br>Crate 2: 81 cm × 81 cm × 91 cm (32 in × 32 in × 36 in);<br>Crate 3: 102 cm × 102 cm × 163 cm (40 in × 40 in × 64 in) |
| Shipping weights for three RGC system crates (approximate)    | Crate 1: 211 kg (464 lb);<br>Crate 2: 123 kg (271 lb);<br>Crate 3: 113 kg (250 lb)  |

<sup>1</sup> Temperatures listed are for models when used with standard pumps; for larger-capacity pump temperatures, contact Lake Shore.

## Ordering information

### RGC Series recirculating coolers

|                |       |
|----------------|-------|
| <b>RGC4-10</b> | 1 W   |
| <b>RGC4-15</b> | 1.5 W |
| <b>RGC4-20</b> | 2 W   |

### Accessories

|                |                                  |
|----------------|----------------------------------|
| <b>10RVP</b>   | Vacuum pumping station           |
| <b>10DDP</b>   | Vacuum pumping station           |
| <b>TS-85-D</b> | Turbomolecular pumping station   |
| <b>336</b>     | Model 336 temperature controller |
| <b>335</b>     | Model 335 temperature controller |
| <b>325</b>     | Model 325 temperature controller |

Copyright © Lake Shore Cryotronics, Inc. All rights reserved.  
Specifications are subject to change.

030724 9:33