


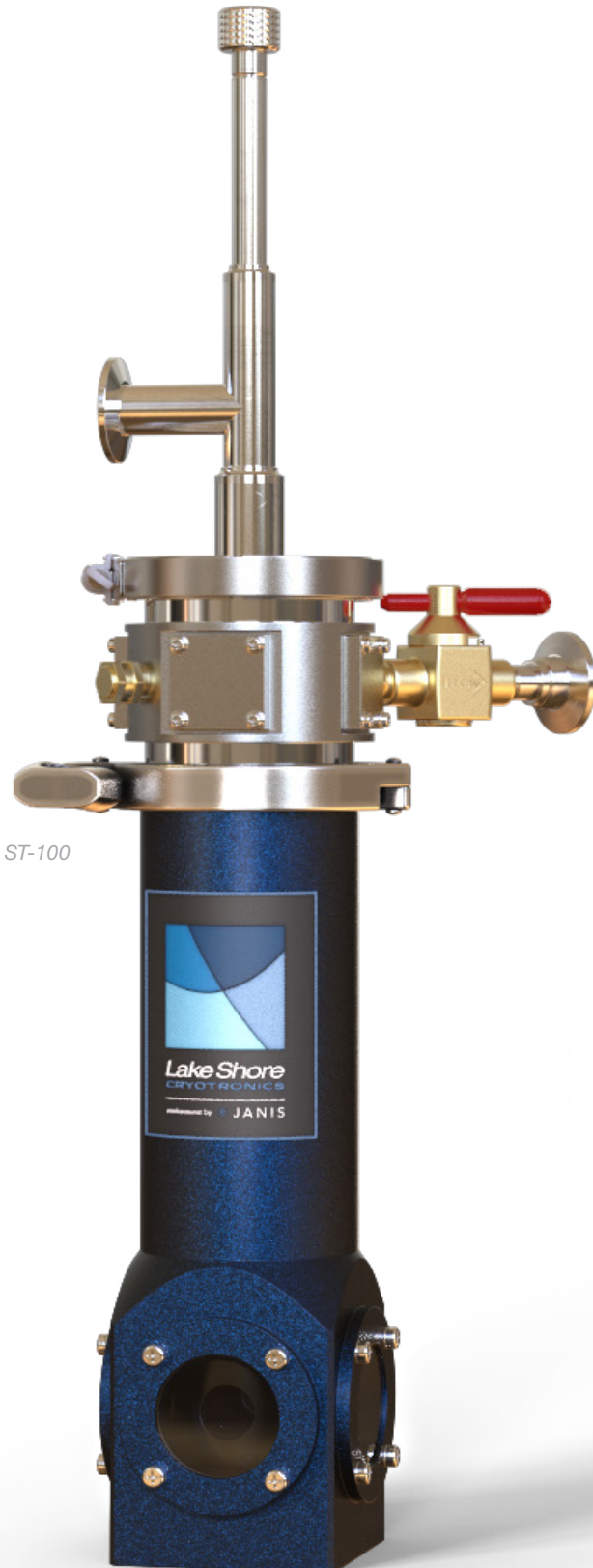

SuperTran Cryostats

ST-100 Series continuous-flow cryostats 2 K to 800 K

ST-100 Series cryostats are continuous-flow cryostats with the sample located in vacuum. They operate with either liquid helium for operation to 2 K or liquid nitrogen for operation to 77 K and withstand high temperatures with options up to 800 K. A high-efficiency transfer line delivers liquid cryogen to the cold finger for cooling, and temperatures below 4.2 K can be achieved with a vacuum pump. By utilizing the built-in heater and temperature controller, these cryostats offer precise variable temperature control within 50 mK.

ST-100 Series cryostats can be combined with Infinite Helium for fully cryogen-free operation throughout the entire temperature range. This enables unattended cryostat operation, making it perfect for extended-duration measurements.

Custom configurations are also available to fit restricted spaces, such as magnet systems or spectrometers, or to accommodate large samples, including semiconductor wafers and “cold plates” that cool multiple samples at once.

Key features

2 K to 500 K (800 K option)

Fast cooldown — 15 min to 5 K

Sample in vacuum

Featured components

High-efficiency, flexible LHe/LN₂ transfer line

Integrated control heater and calibrated control sensor

Polished aluminum thermal radiation shield

ST-100 Series variants

ST-100 maximum temperature = 500 K

ST-100-H maximum temperature = 800 K

Specifications

	ST-100	ST-100-H
Initial cooldown time (LHe to 5 K)	15 min	
Temperature range	2 K to 500 K	2 K to 800 K
Typical temperature stability ¹	±50 mK	
Orientation ²	Any	
Cryogen consumption (LHe room to base temp)	0.4 L	
Cryogen consumption (LHe at 5 K)	0.6 L/h	
Cryogen consumption (LN ₂ at 80 K)	0.1 L/h	
Initial vacuum level requirement ³	~10 ⁻³ Torr	
Typical base pressure during operation	~10 ⁻⁵ Torr	

Size



Height	583 mm (23 in)
Inner diameter (at sample region)	62 mm (2.43 in)
Sample mount diameter	32 mm (1.25 in)
Weight (excluding transfer line)	~4.6 kg (10 lb)
Shipping weight (cryostat only)	8.6 kg (19 lb)
Shipping weight (transfer line)	9.1 kg (20 lb)
Shipping dimensions (cryostat only)	762 × 508 × 508 mm (30 × 20 × 20 in)
Shipping dimensions (transfer line)	2057 × 660 × 127 mm (81 × 26 × 5 in)

¹ Measured with temperature controller

² Cryogen consumption may be higher during non-vertical operation

³ Pressure measured at room temperature prior to adding cryogen

Typical cryogen consumption of optical SuperTran cryostats

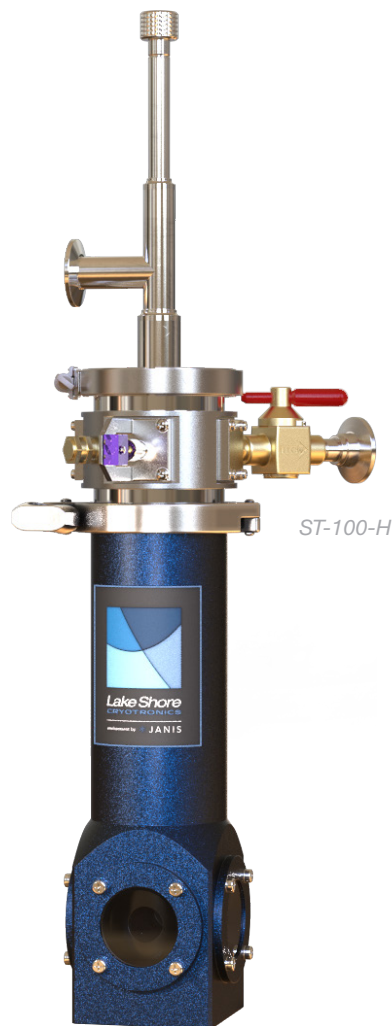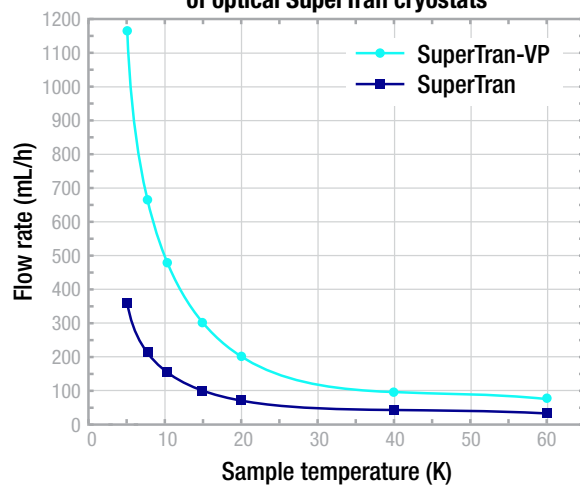

ST-100-H

Complete your setup

Temperature control


Included



Every cryostat includes a Lake Shore temperature controller and calibrated sensor.

MeasureLINK control software

Optional add-on

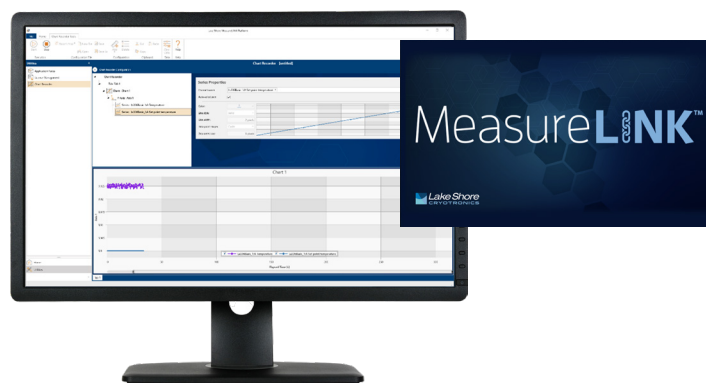

MeasureLINK software enables a wide range of capabilities including charting and logging, system monitoring with a cryostat-specific process view, and controlling Lake Shore equipment as well as third-party instrumentation. No programming required—drag-and-drop to create temperature sweeps, access measurements, and see real-time internal cryostat temperatures in process view.

Source + measure + lock-in


Optional add-on



The Lake Shore M81-SSM provides highly synchronized DC, 100 kHz AC, and mixed DC + AC sourcing and measuring—including both voltage and current lock-in measurement capabilities—for low-temperature material research performed in your cryostat. It supports up to three remote-mountable source and three measure modules per a single M81-SSM-6 instrument and, owing to its modularity, allows signal and source amplifiers to be located as close as possible to the sample being characterized. This minimizes the signal wiring to the sample, reduces noise, and increases measurement sensitivity.

Cryogen-free operation

Optional add-on

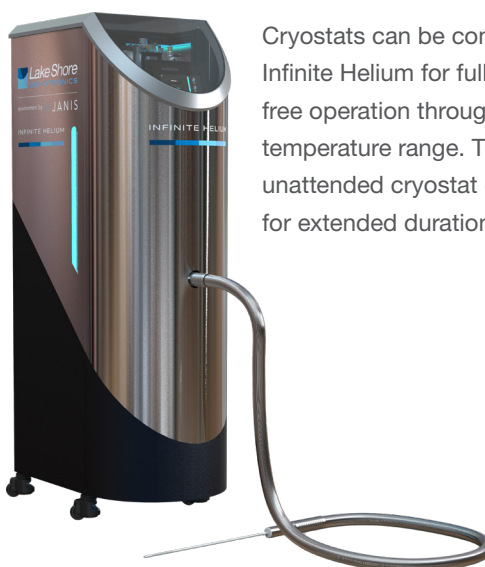

Cryostats can be combined with Infinite Helium for fully cryogen-free operation throughout the entire temperature range. This enables unattended cryostat operation, ideal for extended duration measurements.

Configure your cryostat

1. Select cryostat variant

ST-100	Optical, 2 K to 500 K, calibrated temperature sensor
ST-100-H	Optical, 2 K to 800 K, type E thermocouple
CUSTOM	Custom configurations are available to fit your experiment needs — contact Sales for details

2. Select cryostat configurations

Sample holders

SH-OPTICAL-1.25-STD	Optical
SH-BLANK-1.25-STD	Blank
SH-RESISTIVITY-1.25-STD	Resistivity
SH-FIXED-STD	Fixed
SH-OPTICAL-1.25-800	Optical, high-temperature
SH-BLANK-1.25-800	Blank, high-temperature
SH-RESISTIVITY-1.25-800	Resistivity, high-temperature
SH-FIXED-800	Fixed, high-temperature
CONSULT	DIP

Windows (optical variants only)

See our cryostat window selection guide for additional information.

WR-STD-FS	Fused silica
WR-UV-FS	UV-grade fused silica
WR-STD-SAPH	Sapphire
WR-STD-ZNSE	ZnSe
WR-STD-CAF2	CaF ₂
WR-6MM-KBR	KBr

Mounting flange

BASE-ST-VPF-M	Baseplate for ST-100, ST-300, and VPF-100 Series cryostats — metric threads
BASE-ST-VPF	Baseplate for ST-100, ST-300, and VPF-100 Series cryostats — imperial threads

Optical accessories (ST-100 only)

CONSULT	Reentrant window flange assembly
CONSULT	Bottom window

3. Select pump (optional)

Each cryostat requires a pump to operate. If you do not have an existing pump to use, select one of the pumps below.

10RVP	General-purpose mechanical pumping station
10DDP	General-purpose mechanical pumping station with LN ₂ cold trap and isolation valve
TS-85-D	Turbopumping station

4. Select cryostat wiring

We offer a variety of both unwired and wired feedthroughs to complete your measurement setup. Please refer to the cryostat feedthroughs and wiring guide for more information.

5. Select optional setup configurations

Cryogen-free operation

INFHE-20	Infinite Helium recirculating cooler with base temperature down to <3.3 K
INFHE-15	Infinite Helium recirculating cooler with base temperature down to <3.5 K
RGC4-10	RGC Series recirculating cooler with base temperature down to <4.3 K

Measurement instrumentation

Cryostats come standard with one temperature controller.

M81-SSM electronic synchronous source measure system

Contact us for cables and adapters for M81-SSM/cryostat integration.

M81-SSM-X	M81-SSM instrument with X = 2, 4, or 6 channels; half the channels are dedicated to sourcing and the other to measurement; see modules below
VM-10	AC/DC voltage measure module + lock-in
BCS-10	AC/DC balanced current source module
CM-10	AC/DC current measure module + lock-in
VS-10	AC/DC voltage source module

6. Select optional control software

ML-MCS	MeasureLINK-MCS software with scripting development license; includes lifetime activation for version purchased and full MeasureLINK capability on up to 5 computers with Lake Shore instrument drivers, chart recorder functionality, and drag-and-drop measurement sequences; some application packs sold separately
---------------	--

7. Select additional accessories

Cryostats come standard with one installed temperature sensor. Other sensors are available — contact us.

CX-1050-CU-HT-1.4M	Cernox® magnetic field independent, calibrated
DT-670-CU-HT-1.4H	Silicon diode, calibrated
CONSULT	Thermocouple (ST-100-H only)
CF-100	LHe storage Dewar
LN-50	LN ₂ storage Dewar configured for use with SuperTran cryostats

Copyright © Lake Shore Cryotronics, LLC All rights reserved. Specifications are subject to change.

091126 12:47